Siena



Curved stairlift

Stannah



Make a step change

If a Stannah stairlift could make life easier, why wait?

Life is better when you have the energy to do the things you want to do, like enjoying a visit from the grandchildren or entertaining friends. Our experience has taught us that the sooner a stairlift is installed, the sooner you will be able to start enjoying the things you love doing. Just think how a Stannah stairlift could make your day-to-day life easier.

A stairlift that's right for you and your home.

We know that not all homes are the same. If you think your stairs are unsuitable for a stairlift, we think you might be pleasantly surprised.

- Our stairlifts can accommodate curved staircases of different sizes and layouts.
- We can install our rails on either side of the staircase.
- Our slimline rails and fold-away chairs make the best use of available space, leaving enough room for others to use the staircase.
- All our rails are tailor-made, so they fit your staircase perfectly.
- Contrary to popular belief our rails are fitted to the staircase and not the wall – ensuring minimal mess and bother.
- Installation is quick and easy.
- Our stairlifts can navigate all kinds of corners and landings with a range of rail solutions.

"We couldn't imagine life without our stairlift. It's so easy to keep tidy when one is able to get upstairs to put things away"

Mrs Allsop, Derbyshire



Outside rail



Inside rail



A company you can trust

We are a British, family-run company and have been for five generations. The current managing directors are the great, great grandchildren of Joseph Stannah, who started the company in the 1860s. His engineering principles and strong family values run as true today as they always have.

Now, 140 years later, we have manufactured more than 500,000 stairlifts and installed products in over 40 countries. The half-million milestone was commemorated with a visit from HRH The Prince of Wales – a proud day for all those who work at Stannah.

The Siena

Our Siena stairlift combines fresh, contemporary styling with the best in stairlift engineering. The simple contours along with five classic colour options make the Siena a smart and practical addition to any home.

The Siena comes with a custom-made dual rail that will hug the contours of your stairs precisely and folds away neatly so you can continue to use the stairs as usual if and when you want to.



Standard arm width - 496mm (19.5")



Narrow arm width - 425mm (16.7")

Key features

1. Footrest fold

To store the stairlift away, you can fold away the arms and seat and use the lever to raise the footrest. There's no awkward bending down, which is particularly useful at the top of the stairs. Or for the ultimate ease-of-use opt for our powered footrest button.

2. Swivel seat

Using the levers enables you to swivel the seat easily, making it simpler to get off the chair.

Alternatively, you can choose to have a powered swivel seat, which turns the chair to face the landing automatically.

3. Seatbelt

Choose a seatbelt to suit your style and buckle up for added security (retractable seatbelt shown – optional extra).

4. Safety edges

In the event that something obstructs the stairlift, sensors on the edge of the carriage and footrest will automatically bring the stairlift to a safe stop.



5. Easy-to-use controls

Effortlessly ride the Siena using our simple, intuitive controls that have been designed for those with less dexterity. They work with the lightest of touches, and can be operated using your hand, arm or wrist.





6. Wall controls

Ideal for when you need to summon your stairlift from another floor, or perhaps send it up or down to someone else who wants to use it.




Siena / 8

Using your Siena

The simple-to-use
Siena stairlift has been
designed to make your
life easier and will help
you to rediscover
your independence.
This step-by-step guide
demonstrates some
of its key features.



 The wall-mounted controls at the top and bottom of the stairs make it easy to call the stairlift whether you are upstairs or down.

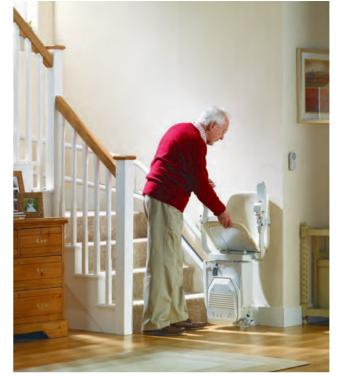

2. Simply unfold the seat and the arms and use the lever to lower the footrest.



3. Once unfolded, make yourself comfortable in the seat. Fasten your seatbelt for extra security.



4. Gently apply pressure to the control in the direction you wish to travel and your Siena will glide smoothly up or down the stairs.



5. When at the top of the stairs, the footrest will stop level with your landing. You can then swivel the chair to face the landing to alight safely.



For added safety, the Siena can remain in the swivelled position to form a barrier at the top of the stairs.



7. When not in use, your Siena folds neatly against the wall, giving you and others plenty of space.

Siena / 10 Siena / 12

Your Siena - tailored for you

Standard Package



In addition to the standard package, you can choose from a range of other optional extras that will ensure your Siena is just right for you. Our sales consultant will guide you through the various options available to you.

Manual swivel seat

Twin-rail pearl colour



Seatbelt



Light-touch control

Left hand Right hand

Upholstery colours and optional extras

Vinyl upholstery












Curved rail colours



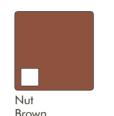










Rail start

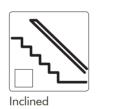







Rail finish





Nicholson Selfe Mobility Services

T: +34 635 726 764 E: info@nsmobilityservices.com



Swivel





Powered

Two-way powered

Seatbelt







Diagonal seatbelt

Five-point harness

Dimensions



Powered footrest button

Min. width of unfolded stairlift: 679mm / 26.7"

Min. width of folded stairlift: 405mm / 16"

Max. weight: 120kg - 135kg / 19 - 21 stone

Standard seat width: 496mm/19.5"

Narrow seat width: 425mm/16.7"